



GETTING STARTED GUIDE

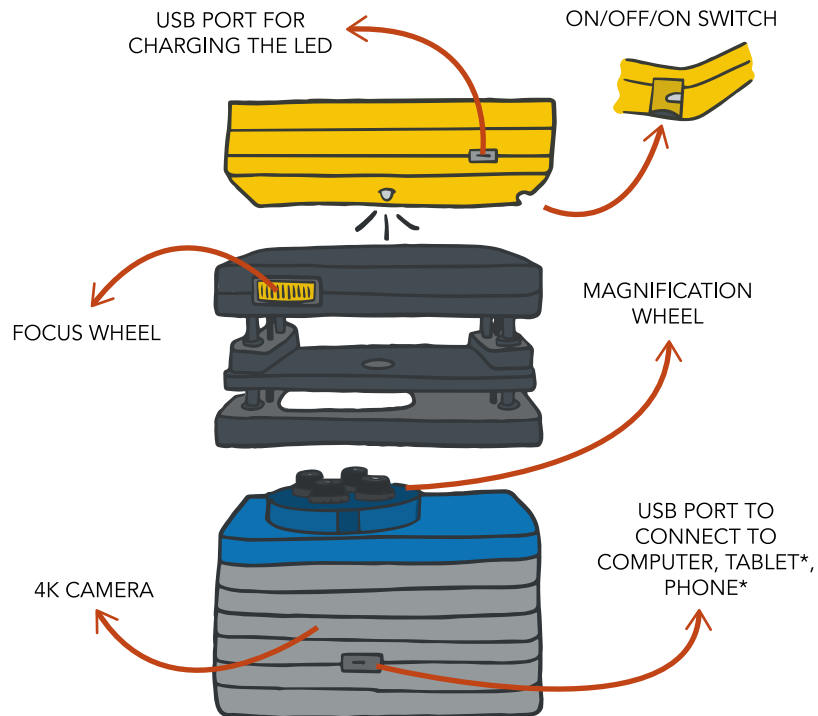

Objectives & magnifications

FIELDS OF OBSERVATION
(width x height)

x2	⇒	4mm x 2.25mm
x5	⇒	2mm x 1.12mm
x10	⇒	1mm x 0.56mm
x20	⇒	0,7mm x 0.39mm

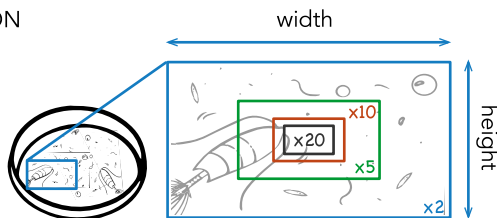

Scan this QR Code:



To access the **resources** for getting started with the Curiosity microscope:

To discover what is being done **all around the world** with the Curiosity microscope:




1. Plankton collection




PLANCTON SCOOP


Dip the **plancton scoop** in the water and draw ∞ to concentrate the plankton

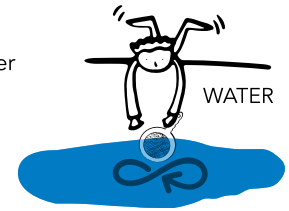

2. Sample preparation

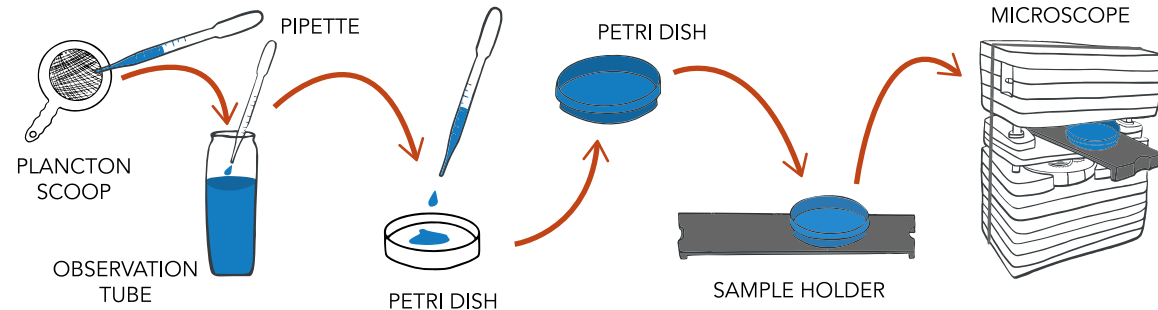

3. Connecting to computer, tablet* or phone*

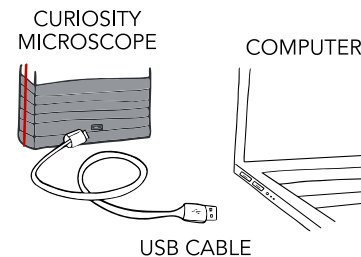

On **Mac**: open the «QuickTime Player» app, select «New video recording», then source «4K USB Camera»



On **Windows**: open the «Camera» app, select the «4K USB Camera» source





On **Android**: install then open the «Next Camera» app, select the «4K USB Camera» source



On **iOS**: install then open the «USB Camera Pro» app, select the «4K USB Camera» source
*only compatible with USBc devices

4. Observations

Turn on the microscope **light** and adjust the sharpness using the **focus wheel**.

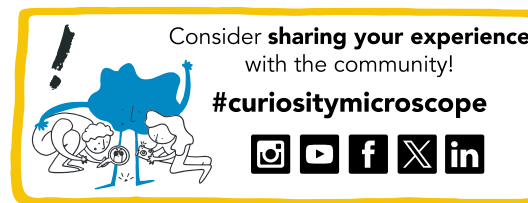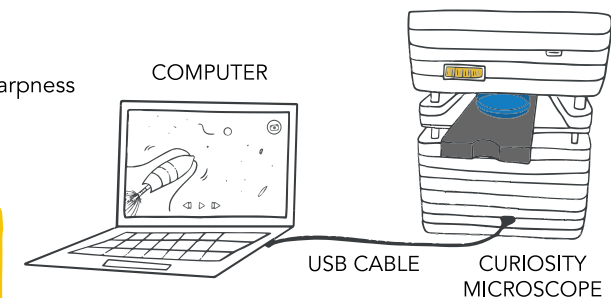

Concept and production SeaLabX – Illustrations Marta Musso